Revive

USER MANUAL | OG357

Hello.

Congratulations on purchasing your new Ogawa Revive OG357.

At Ogawa we understand that every body is different and have carefully designed your massage chair with that in mind. But enough about us, this is about you. Learning how to use your new chair is simple (we designed it that way).


So let's get started.

TABLE OF CONTENTS


 INTRO	Dimensions	2
	Contents	3
 CHAIR	Move & Clearance	4
	Components	5
	Massage Points	6
	Specification	7
 MESSAGE	Operation Guide	8-10
	Pre-Programed Description	11-12
 CARE	Safety	13-14
	FCC Statement	15
	Grounding Instruction	15
	Troubleshooting	16
	Care and Storage	16
	Warranty	17

DIMENSIONS


First, let's get your chair inside.

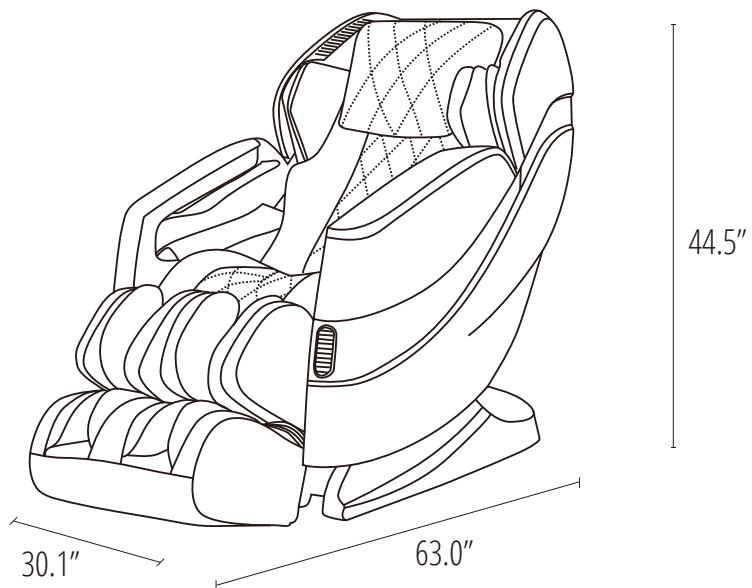Box 1

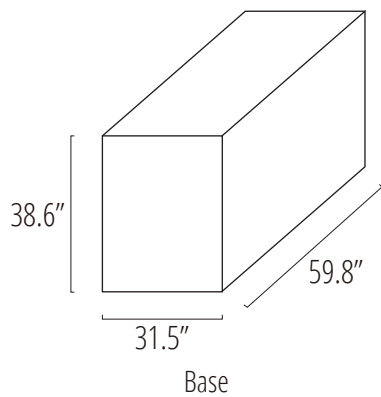Box 2

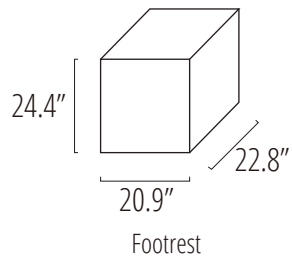CONTENTS

Everything you need for a relaxing experience.


User Manual


Installation Manual

Warranty Postcard

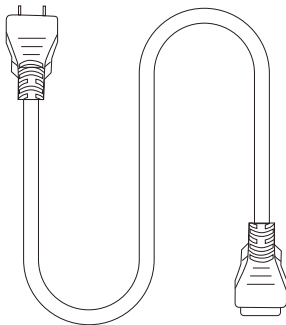Power Cord

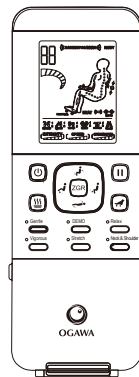Remote Control


OG357 Massage Chair


Move & Clearance

A few simple steps to find the best place for your new chair.

HOW TO MOVE

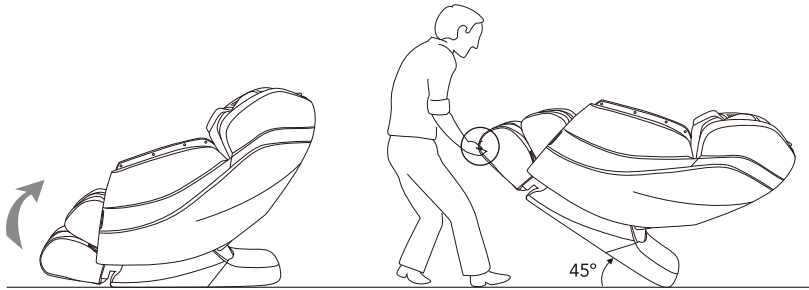
Lift up footrest 45 degrees onto the wheels

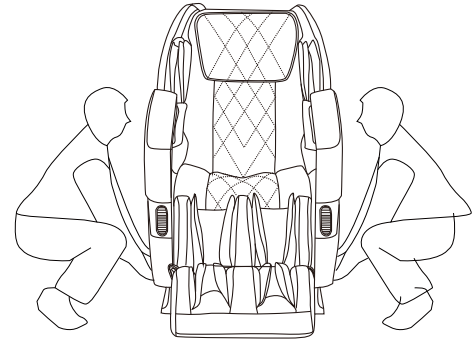
Two people are needed to lift the chair

CLEARANCE

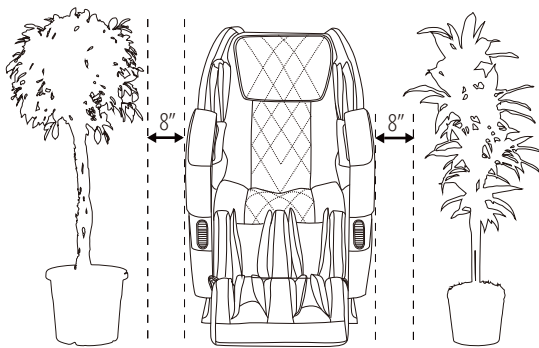
Reserve space on left and right

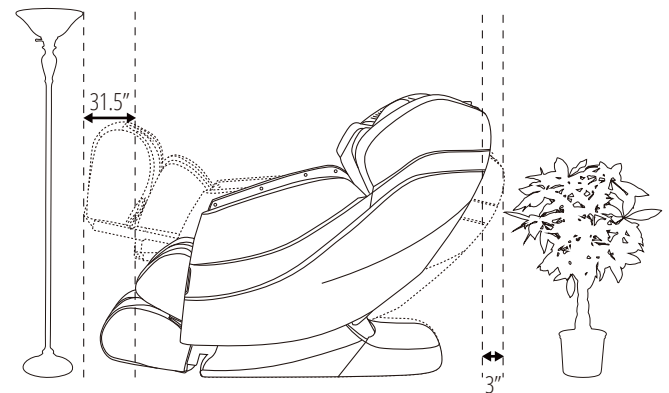Reserve space on the front and back

CHAIR COMPONENTS

Some of your chair's most important features.

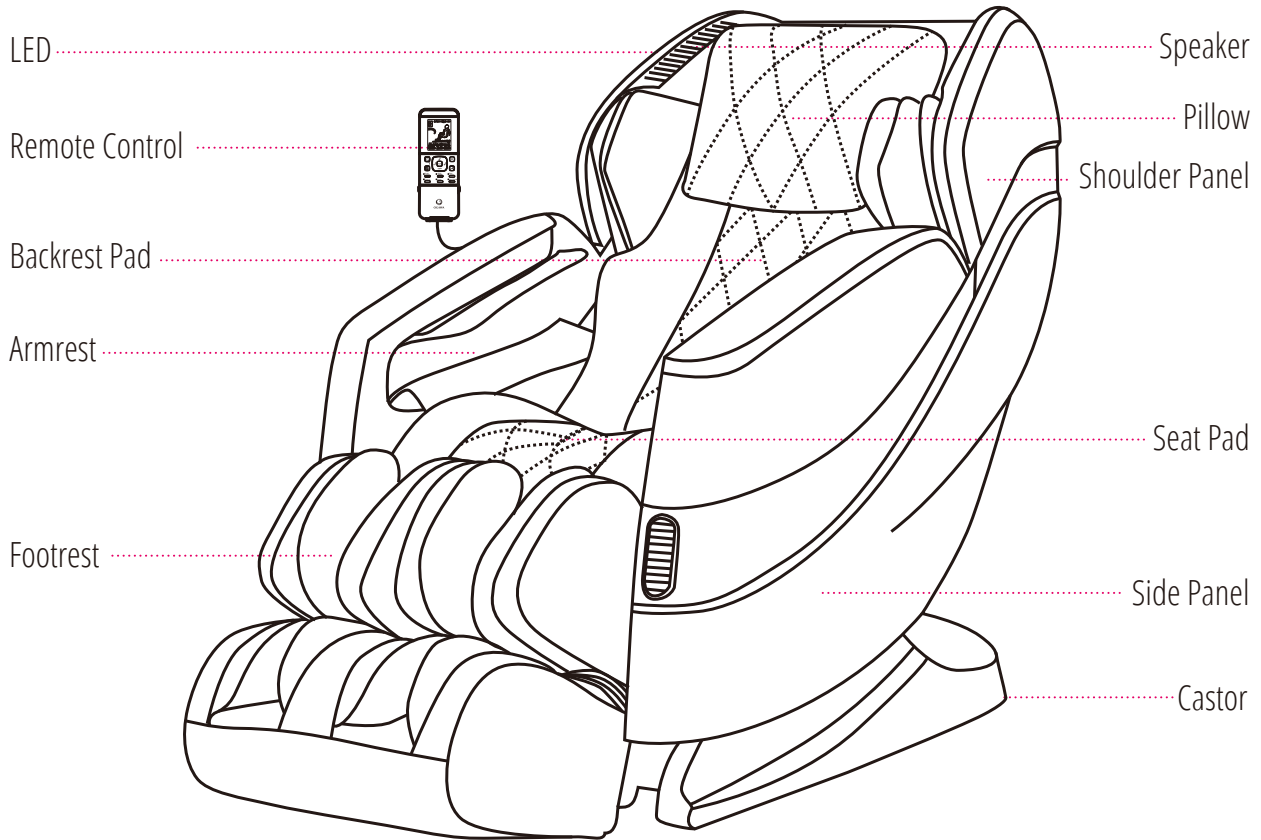
Note: The speaker will only operate when the massage chair is in operation

To play music using your massage chair's speakers please connect your smartphone, tablet, or other electronic device to your chair using bluetooth.

MASSAGE POINTS

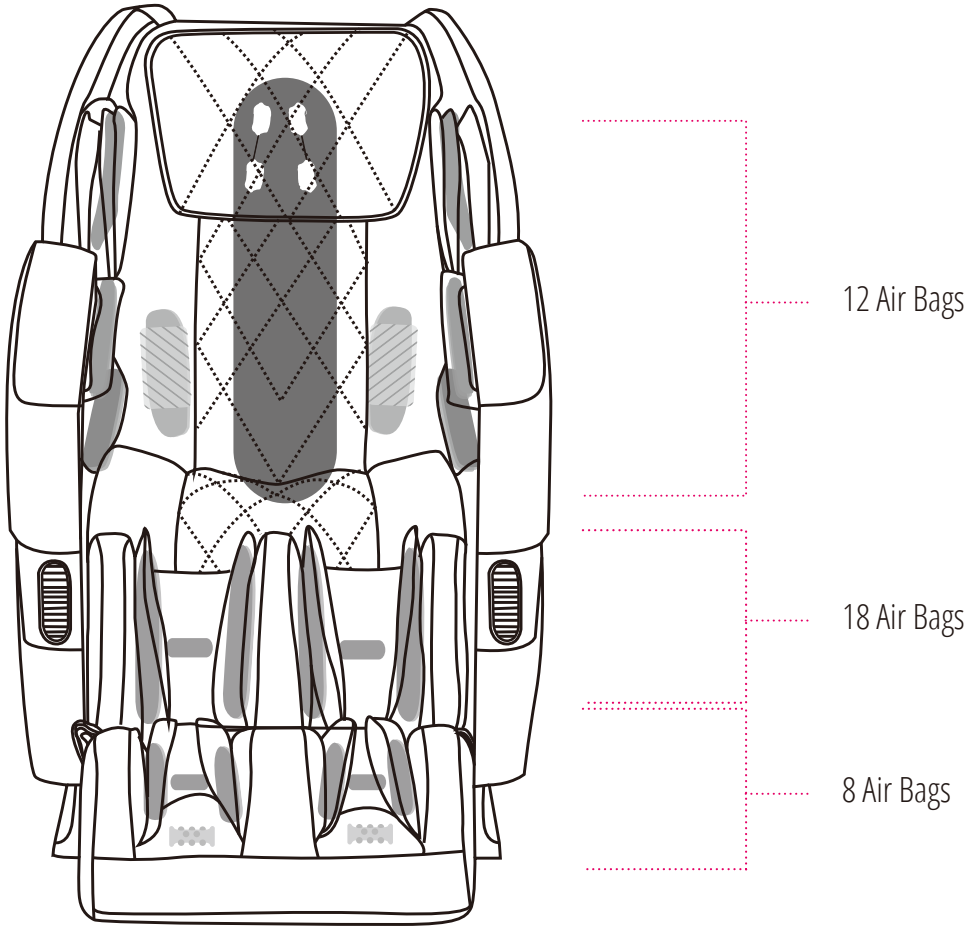- Heaters
- Back Massage Mode
- Air Bag
- Roller

CHAIR SPECIFICATIONS

A few important facts about your massage chair.


Model No.

OG357

Max Massage Time

30 minutes

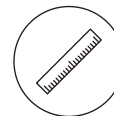Dimensions (L x W x H)

Chair dimensions: 63.0"x30.1"x44.5"
Reclined dimensions: 72.1"x30.1"x37.0"


Shipping Weight

Gross weight: 225 lbs
Net weight: 185.2 lbs

Usage Condition

Environment temperature: 10°C-40°C
Contrasting humidity: 30-85RH

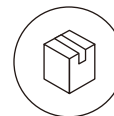Storage Condition

Storage temperature: -20°C-60°C
Storage humidity level: 30-85RH

Power Consumption

100W

Rated Voltage

110 - 120V

Rated Frequency

60Hz

Power of Speakers

14W (2 sides)

Length of Wiring

Controller wiring 27.6"
Power supply wiring 70.9"

Safety Feature

Equipped with overheating and power surge detectors safety feature

OPERATION GUIDE

Connect to the power supply, and turn the switch to “I” position, chair is now in idle mode.

During its idle mode ON/OFF button will blink.

During its idle mode, only ON/OFF, Pause, Leg rest Recline/Restore, Chair Recline/Restore buttons can be operated.

1. ON/OFF

The Start/Stop Control Button

Press this button to turn on the massage chair. Press any massage function button to start a massage. Press this button again to stop all massage functions and turn off the chair.

2. PAUSE

Press this button once to stop current massage function and countdown timer. During the pause mode, the pause button will be blinking. Press this button again to cease the "Pause" function and resume previous massage function, as well as countdown timer.

3. HEAT

Lumbar Heat function Start/Stop Control Button

Press this button to activate lumbar heat function, press this button again to stop the heat function.

Note: During any auto program, pressing this button can also disable/enable the heat function.

4. FOOT ROLLING

Foot Rolling Function Start/Stop Control Button

Press this button once to activate rhythmic rolling massage at soles, press this button again to stop rolling massage at soles.

Note: During an auto program, pressing this button can also disable/enable the foot rolling function.

5. GENTLE

A program designed to provide a quiet, gentle massage that is good for beginners or those that prefer a soft touch.

6. DEMO

This message will quickly demonstrate all of the massage features in 5 minutes.

7. RELAX

A program designed to provide a mild massage that is great for deep relaxation, best used before going to sleep.

8. VIGOROUS

This message will combine vigorous tapping on your back and an air massage focused on your arms and legs.

9. STRETCH

This message will stretch your legs and back with a combination of gentle tapping, kneading, rolling, and air compression.

10. NECK & SHOULDER

This message will focus on your shoulders and neck using a combination of keading and shiatsu techniques.

11/12. LEGREST RECLINE AND RESTORE

Angle Adjustment Button for Leg Rest

Press and hold either button until the leg rest reaches your desired position.

13/14. CHAIR RECLINE AND RESTORE

Hold this button down to comfortably lay flat or to restore the chair to its original position.

15. ZGR

Press this button in order to active the zero gravity position, where your feet are raised above your heart.

16. LED LIGHTING

LED

LED Lighting Function Start/Stop Control Button

Every press of this button will change the colors of LED lighting in the sequence of "Red > Blue > Green > Three color auto cycle > Close".

Note: When the massage chair is powered on, the three colors auto cycle, and the auto programs have preset LED colors. You can also adjust the LED color by pressing the LED button.

17. SPOT/PARTIAL

Spot or Partial Back Massage

Press this button once to activate localized spot massage; the original massage mode will continue to run. Every press of this button will change the massage function in the sequence of "SPOT > PARTIAL > STOP" in cycle.

Note: Partial Massage will start at the spot chosen and will move up and down 8 cm. Spot and partial are only effective in manual massage mode. Spot function will not work while rolling function is activated.

18/19. UPWARD/DOWNWARD ADJUSTMENT

When the spot and partial massage function are activated, this button can be used to move the massage rollers. Hold this button to move the massage rollers upward/downward respectively. Release the buttons when reaching desired position.

Note: When the auto program is engaged, it will scan the back length automatically which will set the highest position. The position can be adjusted by pressing these two buttons after hearing beep.

20. MANUAL MASSAGE FUNCTION

M

KNEADING, SWEDISH, TAPPING, SHIATSU, CLAPPING, ROLLING

Press this button to activate the kneading manual massage function for a full back massage at medium strength. Every time this button is pressed, the manual massage function will change in the sequence of Kneading > Swedish > Tapping > Shiatsu > Clapping > Rolling > Stop.

Note: Kneading/Swedish/Tapping/Shiatsu/Clapping have two modes.

21. STRENGTH

The Strength Adjustment Button for Kneading Function

When the kneading function is engaged, every press of this button will change the massage speed strength interval in the sequence of "Lv1 > Lv2 > Lv3 > Lv4 > Lv5".

22. SPEED

The Speed Adjustment Button for Tapping, Clapping, Swedish and Shiatsu Function

When Tapping, Clapping, Swedish or Shiatsu function is engaged, every press on this button will change the massage speed interval in the sequence of "Lv1 > Lv2 > Lv3 > Lv4 > Lv5" in cycle.

23. WIDTH

The Width Adjustment Button for Tapping, Clapping, Shiatsu and Rolling Function

When Tapping, Clapping, Shiatsu or Rolling function is engaged, every press of this button will change the massage width interval in the sequence of "Narrow > Medium > Wide" in cycle.

24. VIBRATION

Seat Vibration Function Start/Stop Control Button

Press this button once to activate rhythmic seat vibration, press this button again to stop the vibration function.

Note: During an auto program, pressing this button can disable/enable the vibration function.

25. FULL AIR

Press this button to activate the air massage function for the full body (Calf, Foot, Arm, Shoulder). Press this button again to stop the function.

26. PARTIAL AIR

Press this button to activate air compression massage from shoulder to soles. Every time this button is pressed, partial air massage function will change in the sequence of Calf & Foot > Arm > Shoulder > Stop.

Note: If you press the full air button, it will switch from manual air function to full air function.

27. INTENSITY

The air compression intensity adjustment button has 5 levels of intensity. When the air massage function is engaged, every press of this button will change the air intensity in the sequence of "Lv1 > Lv2 > Lv3 > Lv4 > Lv5" in cycle.


PRE-PROGRAMED DESCRIPTION

GENTLE

A program designed to provide a quiet, gentle massage that is good for beginners or those that prefer a soft touch.


FOCUS AREA

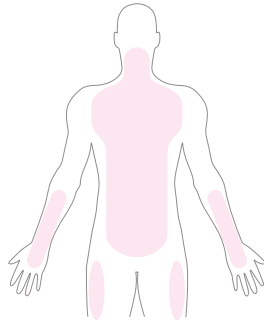INTENSITY

FOCUS TECHNIQUE

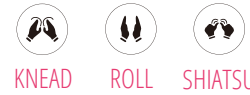
DEMO

This massage will quickly demonstrate all of the massage features in 5 minutes.

FOCUS AREA

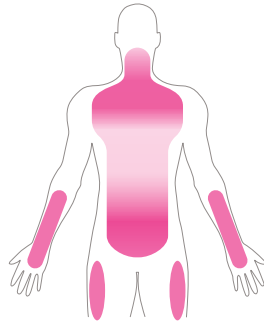INTENSITY

FOCUS TECHNIQUE


STRETCH

This massage will stretch your legs and back with a combination of gentle tapping, kneading, rolling, and air compression.

FOCUS AREA

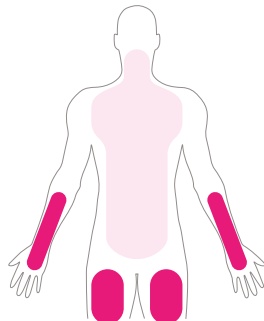INTENSITY

FOCUS TECHNIQUE


RELAX

A program designed to provide a mild massage that is great for deep relaxation, best used before going to sleep.

FOCUS AREA

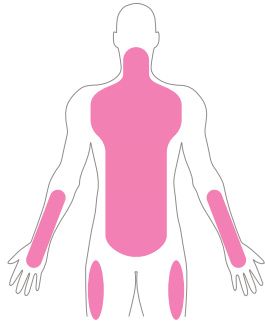INTENSITY

FOCUS TECHNIQUE

KNEAD

ROLL


CLAP

Neck & Shoulder

This massage will focus on your shoulders and neck using a combination of kneading and shiatsu techniques.

FOCUS AREA

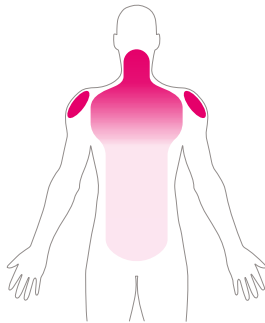INTENSITY

FOCUS TECHNIQUE

KNEAD

ROLL


SHIATSU

Vigorous

This massage will combine vigorous tapping on your back and an air massage focused on your arms and legs.

FOCUS AREA

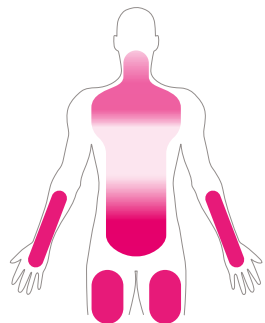INTENSITY

FOCUS TECHNIQUE

KNEAD

CLAP

TAP

SAFETY PRECAUTIONS

We're looking out for you.

IMPORTANT SAFETY INSTRUCTIONS

When using the massage chair, basic precautions should always be followed, including the following:
Read all instructions before using this massage chair.

DANGER—To reduce the risk of electric shock:

- Always unplug this massage chair from the electrical outlet immediately after using and before cleaning.
- Never use pins or other metallic fasteners with this appliance.
- Carefully examine the covering before each use. Call for service if the covering shows any sign of deterioration, such as blistering, peeling or cracking.
- Keep dry—Do not operate in a wet or moist condition.

WARNING—To reduce the risks of burns, fire, electric shock or injury to persons:

- Read all instructions before using.
- A massage chair should never be left unattended when plugged in. Unplug from the outlet when it is not in use and before putting on or taking off parts.
- Close supervision is necessary when this massage chair is used by or near children or disabled persons. Keep children away from reclining backrest and extended leg rest.
- Use this massage chair only for its intent use as described in this manual.
- Never operate this massage chair if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Contact Customer Service Department.
- Keep the cord away from heated surfaces.
- Never operate the massage chair with the air openings blocked. Keep the air openings free of lint, hair and the like.
- Never drop or insert any object into any opening.
- Do not use outdoor.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- To disconnect, turn all controls to the off position, then remove the plug from the outlet.
- Do not stand on or in the massage chair. Use only while seated.
- The product is not intended for use by children or persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the product by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the product.
- Temperatures sufficiently high to cause burns may occur regardless of the control setting. Do not use on an infant or invalid or on a sleeping or unconscious person. Check the skin in contact with the heated area of the appliance frequently to reduce the risk of blistering.
- Burns may result from improper use.
- Do not crush in order to avoid sharp folds.
- Keep children away from extended foot support (or other similar parts).
- DO NOT WET.
- DO NOT USE PINS.
- NEVER REMOVE COVER.
- Connect this appliance to a properly grounded outlet only. See Grounding Instructions.

SAVE THESE INSTRUCTIONS

This massage chair is intended for household use only.

1. Usage Environment

- Do not use the massage chair in an excessively humid or dusty environment as this may result in malfunction or electric shock.
- Do not use the massage chair in a room with a temperature of 104°F or higher.
- Do not expose the massage chair to heaters, stoves or direct sunlight.
- Place and use the massage chair on an even, non-slip floor.
- The massage chair is designed for home use and not for commercial purposes.

Caution

- This massage chair is to be used indoors only.
- If you use the massage chair in a cold room, do not increase the room temperature abruptly. It is recommended that you increase the temperature gradually to a normal level.
- If the massage chair is stored in a cold place and is being brought into a warm environment, it is advisable to wait for one hour before use. This is because its performance may not be at optimum level, as there may be water droplets on the mechanical parts due to consideration. Using the massage chair under such conditions may result in malfunction.

2. Safety Precautions

- To ensure safe and correct use of the massage chair, do not operate it without reading this operation manual carefully.
- Do not use the massage chair when the upholstery or seat is removed. This may cause injury or malfunction.
- Do not use the massage chair in combination with other therapeutic equipment or electric blanket etc., as this may result in ineffectiveness or injury.
- Do not use the massage chair an hour before or after eating as this may have adverse effects on the user.
- Do not allow children or pets to play around the massage chair, e.g. behind the backrest, under the seat base or legrest of the massage chair. This may result in injury.
- Do not use the massage chair when your body is wet.
- Do not operate the massage chair with wet hands.
- It is not recommended to use the massage chair for more than 30 minutes.
- Do not perform continuous massage on the same spot of your body for more than 5 minutes at a time as this may result in excessive stimulation and may have adverse effects.
- Do not insert your hand or foot along the paths of the massage rollers during use as this may cause injury.
- Stop operation immediately if you feel any discomfort.
- This product is not intended for self-treatment of conditions that should be managed by a qualified health care provider.
- People on medication or with medical conditions, please consult your doctor before use. Please do not use the product if you are not feeling well.

3. People with Medical Conditions

Consult your doctor before use if you:

- are under medical rest as ordered by a doctor.
- have spinal disorders, an abnormal spinal condition or have suffered a spinal injury.
- have back problems.
- have diabetes, osteoporosis or sensory impairment.
- have joint dysfunction such as rheumatism, hammertoe or gout.
- have a pacemaker or other electronic medical devices.
- may be pregnant.
- have phlebitis or thrombosis.
- have an increased risk for blood clots.
- had recent surgery.
- Have surgical pins, screws, or anything mechanical in the legs, ankles or feet.

Caution

- This massage chair is not intended for use by children or persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction or instruction concerning use of the massage chair by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the massage chair.
- Avoid using the massage chair under the influence of alcohol.
- Do not apply massage rollers directly on the head, elbow or knee joints, torso or abdomen.
- If your legs or feet slip off the leg rest during massage, do not force them back in as this may result in injury.
- Remove the power plug from the electrical outlet immediately should liquid be accidentally spilled onto the massage chair. Do not use wet hands to connect or disconnect the power plug from the electrical outlet.
- Do not tamper with the power cord or use the massage with a damaged cord as this may result in malfunction and electric shock.

4. Instances When the Massage Chair Should Not Be Used

In the following events, please switch off the main power switch immediately and disconnect the power cord from the electrical outlet:

- If liquid is accidentally spilled onto the massage chair. This may result in electric shock or malfunction.
- If the fabric of the backrest is torn and the internal components are exposed.
- If pain or discomfort is felt during massage, stop operation immediately and consult your doctor.
- If you detect any malfunction or any other abnormal condition during operation.
- If there is power failure. Injury may occur when power is restored unexpectedly.
- If there is lightning.

5. Assembly and Repair of the Massage Chair

- Do not disassemble the backrest cover of the massage chair. Touching the internal components may result in malfunction or electrical shock.
- This massage chair must only be professionally serviced or repaired. Do not attempt to disassemble or repair the massage chair by yourself.

6. Things to Be Noted About the Power Plug and Cord


- Check that the Alternating Current (AC) voltage corresponds with the specifications indicated for the massage chair.
- Do not connect or disconnect the power plug from the electrical outlet with wet hands. This may result in malfunction or electric shock.
- When disconnecting the power plug, pull it by holding the plug, not the cord.
- Although operation will be stopped automatically by the Auto-Timer function, always remember to switch off the main power after use.
- Do not place the power cord under the massage chair or any other heavy objects.
- Do not wind the power cord around the massage chair as this may damage the cord and result in a fire or electric shock.
- Do not operate the massage with a damaged cord or extension cord.
- Do not use the massage if the electrical outlet is loose.
- If the power cord is damaged, it should be replaced by a qualified technician.

GROUNDING INSTRUCTIONS

This product must be grounded. If it should malfunction or break down, grounding provides a path of least resistance for the electric current to reduce the risk of electric shock. This product is equipped with a cord with an equipment-grounding conductor and a grounding plug. The plug must be plugged into an appropriate outlet that is properly installed and grounded in accordance with all local codes and ordinances.

Danger-Improper connection of the equipment-grounding conductor can result in a risk of electric shock. Check with a qualified electrician or service technician if you provided with the product-if it will not fit the outlet, have a proper outlet installed by a qualified electrician.

This product is for use on a nominal 120V circuit and has grounding plug that looks like that plug illustrated below. Make sure that the product is connected to an outlet having the same configuration as the plug. No adapter should be used with this product.

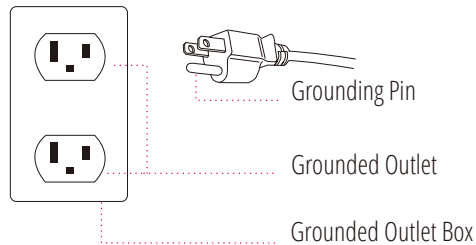TROUBLESHOOTING

This massage chair is controlled by a microcomputer and is thus susceptible to power surges of +10 to-15%. The following situations may arise in the event of a power surge:

Situation	Solution
The massage chair is connected to the power supply but is not functioning	<ol style="list-style-type: none">1. Check to see that the plug is properly connected and working2. Ensure the massage chair switch is in the 'I' position
The remote control works but the massage chair is not functioning.	Possible overheating due to prolonged use. Turn off the power for at least 30 to 50 minutes
Both the remote control and the massage chair are not functioning	<ol style="list-style-type: none">1. Ensure the chair is plugged into a working outlet2. Turn the power OFF and ON again

If you continue to have problems, please read the warranty card and contact the Cozia service department.

PRODUCT CARE AND STORAGE

1. Storing the Massage Chair

- Wipe off dust and dirt from the massage chair with a dry, soft, clean cloth after use.
- Keep the massage chair clean and away from places of high humidity.
- If the massage is not to be used for a long period of time, cover the massage chair with a dust cover and unplug its power plug from the electrical outlet.

Caution

- Do not store the massage chair in direct sunlight or in places of high temperature. This may cause the color of the upholstery to fade.
- Do not allow the massage chair to remain in contact with vinyl or wallpaper for extended periods of time as this could lead to discoloration.

2. Cleaning the Massage Chair

- Before cleaning, ensure that the main power switch is off and the power plug is disconnected from the electrical outlet.
- Wipe off dust and dirt from the massage with a dry, soft, clean cloth after each use.
- This product can be cleaned thoroughly with a slightly damp cloth.
- Do not iron the upholstery.
- Do not wash the upholstery using a home washing machine as this may result in deterioration of the upholstery's quality.

Caution

- Do not use abrasive cleaning agents like benzene or thinner to clean the massage chair. This may cause the upholstery to fade.

3. Cleaning the Controller

- Use only a clean, dry cloth to clean the controller.

Caution

- Do not use a wet cloth to clean the controller. This may result in malfunction.

COZZIA CONSUMER WARRANTY

Ogawa products have the following warranty:

This Ogawa massage chair provides a one year in-home service warranty and a two-year parts warranty. Ogawa customer service will not issue Return Material Authorizations (RMAs) for products. Ogawa products will be either repaired by the consumer or by an in-home technician. Ogawa will not issue Return Material Authorizations (RMAs) for buyer's remorse.

Service And Technical Support

Consumers may contact Ogawa's customer service department at 844-438-6429 between 9:00 am and 5:00 pm PST, Monday through Friday or email us at support@ogawaworldusa.com or visit our web at www.Ogawaworldusa.com for warranty or service issues.

Consumers are required to provide the unit serial number and dated proof of purchase (sales receipt) when they contact the Ogawa service department regarding a repair.

A Ogawa customer service representative will attend to most consumer inquiries, but in some cases a technical service specialist will provide advanced support.

Non-warranty repair is provided on a "per incident" basis. Ogawa customer service will verify that the unit has failed and provide instructions for repairing a unit. All applicable repairs, parts, shipping, handling, local tax and a "per incident" fee will be charged for non-warranty repairs and support calls.

Proof of purchase (original receipt) is required for all warranty repairs or service.

Extended Warranties:

Ogawa offers extended warranty programs to enhance ownership of Ogawa product. Extended warranty provides: 3-year limited repair or replacement massage chair warranty with on-site service. Please visit www.ogawaworldusa.com for the latest terms and pricing for extended warranties.

Product Repairs

Many Ogawa product repairs may be performed by the consumer in their own home with parts and instructions supplied by Ogawa such as remote controller replacements that don't require a technician. Canadian consumers may contact Ogawa's customer service department at 844-438-6429 between 9:00 am and 5:00 pm PST, Monday through Friday for warranty or service issues.

- Floor Models and Demonstration Units have a prorated warranty period that begins the day the model is placed on the floor and connected to an electrical outlet by the retail partner. Proof of purchase is required in order to obtain warranty service and parts. Non-warranty parts and service are available at a "per incident" charge.
- Warranty Coverage Period - Warranty coverage starts at the time that the consumer PURCHASES their chair. Proof of purchase is required for all warranty repairs. The warranty time period is measured by continuous calendar days based on a seven-day week.
- Field Service – Ogawa customer service department will diagnose and schedule an authorized service provider to repair the chair in the consumer's residence.
- Parts - All parts deemed necessary for repairs will be shipped to the consumer at no charge by the Ogawa service department.

CONSUMER WARRANTY LIMITATIONS AND EXCLUSIONS

- Warranties are only valid in the 48 contiguous United States and Canada and are contingent upon the consumer operating the chair according to the corresponding instruction manual.
 - The warranty only applies to Ogawa products and does not include any accessories or enhancements.
 - Field service requires pre-approval and must be performed by Ogawa's authorized field service personnel. Field service is only available in the United States and in major metropolitan areas in Canada.
 - Softening of foams and filling composites in pillows, pads, and memory foam occurs naturally and does not constitute a defect under this warranty.
 - Fading, wear and piling of fabrics occurs naturally and does not constitute a defect under this warranty.
 - Natural markings, grain and dye variations in leather do not constitute a defect under this warranty.
 - Warranties do not cover any loss or damage resulting from improper installation, unauthorized repairs or modifications, improper use of electrical/power supply, loss of power, dropped product, a malfunction or damage of an operating part from failure to provide manufacturer's recommended maintenance, transportation damage, theft, abuse, misuse, neglect, vandalism or environmental conditions (fire, floods, rust, corrosion, sand, dirt, windstorm, hail, earthquake or exposure to weather conditions), loss of use during the period the product is at a repair facility, or otherwise awaiting parts or repair.
 - Damage due to shipping and handling does not constitute a defect under this warranty.
 - Under no circumstance shall Ogawa or its representatives be liable for indirect, consequential, or incidental damages (including damages for lost profits, business interruption, bodily injury, medical and the like), even if any party has been advised of the possibility of such damages.
 - Warranties are non-transferable and shall be in lieu of any other warranty, express or implied, including but not limited to any implied warranty or merchantability or fitness for particular use.
 - Ogawa's sole liability and the purchaser's exclusive remedy shall be for the repair, or at Ogawa's option, for the replacement of the defective part.
Notwithstanding the above, if replacements parts for defective materials are not available, Ogawa reserves the right to make substitutions in lieu of repair or replacement.
-
- All warranties begin on the DATE OF PURCHASE; no allowance or extension is offered for delivery and/or installation.
 - Warranties do not apply to rental, business, commercial, institutional or other non-residential users.

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